

S **semper**
med

supreme duo

**2 PAIRS OF GLOVES
ONE POUCH**

ONE STOP DOUBLE DONNING SOLUTION

This complete double donning solution comes double-packed with one pair of sempermed® supreme** and one pair of sempermed® supreme green. The perfect supreme couple for enhanced safety and fast and clear detection of perforations during surgeries.





supreme duo




Fast and clear identification of punctures

Our optimized transparent overglove makes punctures more visible. The green underglove, combined with the transparent overglove, clearly shows any fluid penetration between the layers.



Optimized sensitivity

The reduced wall thickness of the sempermed® supreme green underglove, used with the sempermed® supreme overglove, perfectly balances sensitivity and surgical finesse with supreme safety.

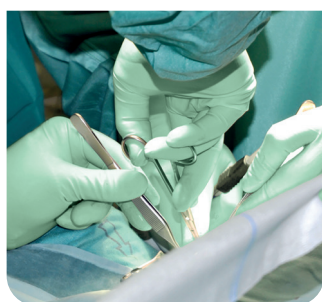

Easy to don

The unique sempermed® quick-don inner coating of both sempermed® supreme gloves makes them easy to don and change during surgeries.



Comfort & fit

Utilizing ergonomically designed formers that provide hand-specific fit for both sempermed® supreme gloves reduces hand fatigue and stress, offering more comfort and a better fit during surgeries.



Technical data

Type	sterile surgical glove for single use, powder free, with synthetic inner lining
Color	natural white / green
Size / Overall length as per EN 455-2	5½, 6 and 6½ / median 270 mm 7, 7½ and 8 / median 280 mm 8½ and 9 / median 285 mm
Wall thickness measured in single layer in the palm area	0.19 - 0.23 mm / 0.17 - 0.20 mm
Barrier performance	AQL 0.65
Material	natural rubber latex (NR)
Glove shape	fully anatomical with rolled rim
Force at break as per EN 455-2	≥ 12 N, exceeding EN 455-2 (≥ 12 N instead of 9 N)
Shelf life in original packaging if stored according to product specification	3 years
Sterilization	radiation STERILE R
Labeling & conformity to standards	EN 455-1/-2/-3/-4, EN ISO 21420, EN ISO 374-1 (Type B), EN ISO 374-2, EN 16523-1, EN ISO 374-4, ISO 374-5, EN 421, EN 556, ISO 2859-1, ISO 15223-1, ASTM F1671/F1671M

Purpose

Single-use medical device class IIa
according to MD Directive 93/42/EEC



Single-use protective glove – PPE category III
(protection against chemical substances for limited time)*
according to PPE Regulation (EU) 2016/425



Key benefits

Fast and clear identification of punctures
Optimized sensitivity – thinner wall thickness of underglove
Unique sempermed® inner coating – easy to don and change during surgeries

Packaging

- Left and right glove turned inside out, in inner pouch, ozone-tight, sealed in medical fiber-free peel pack, double-packed.
- In dispenser carton with sterilization indicator: 25 pouches
- In transport carton with sterilization indicator: 150 pouches

Article numbers

size 5½	3000007261	size 7½	3000007265
size 6	3000007262	size 8	3000007266
size 6½	3000007263	size 8½	3000007267
size 7	3000007264	size 9	3000007268